EN

User Manual

for laypersons and professional users

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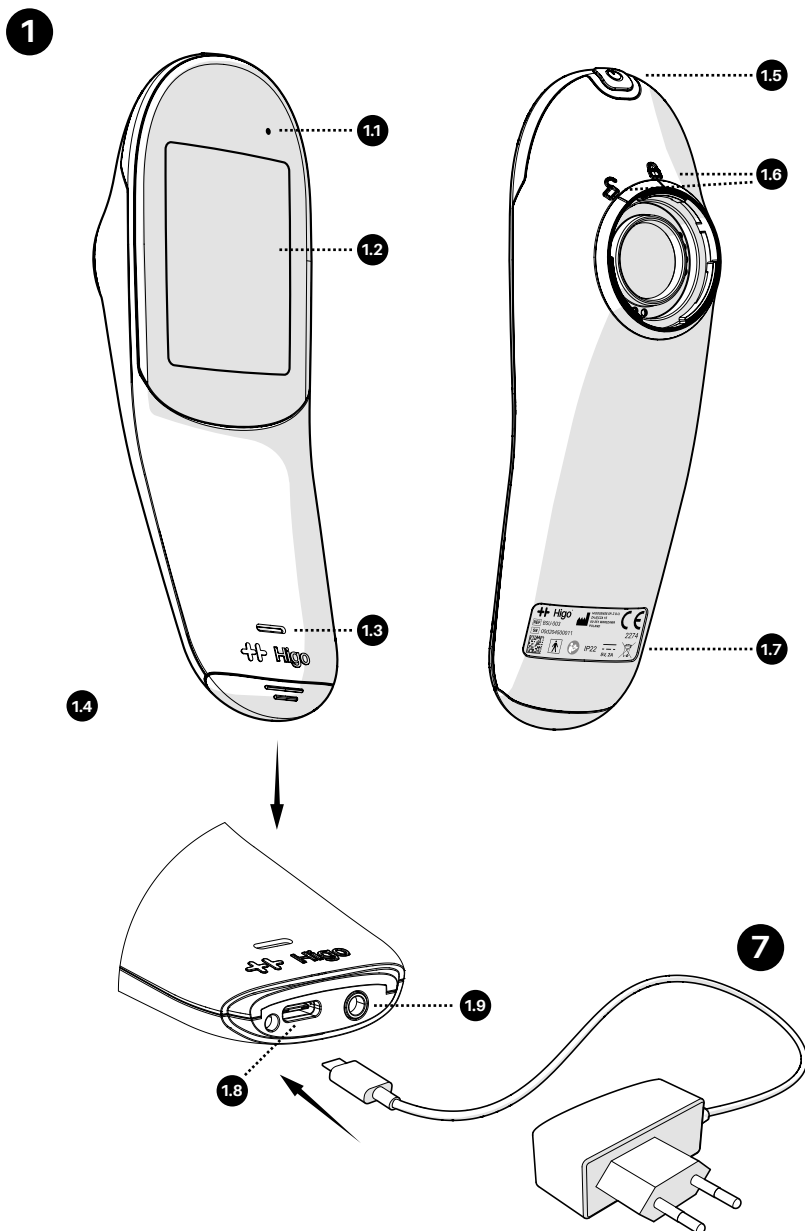
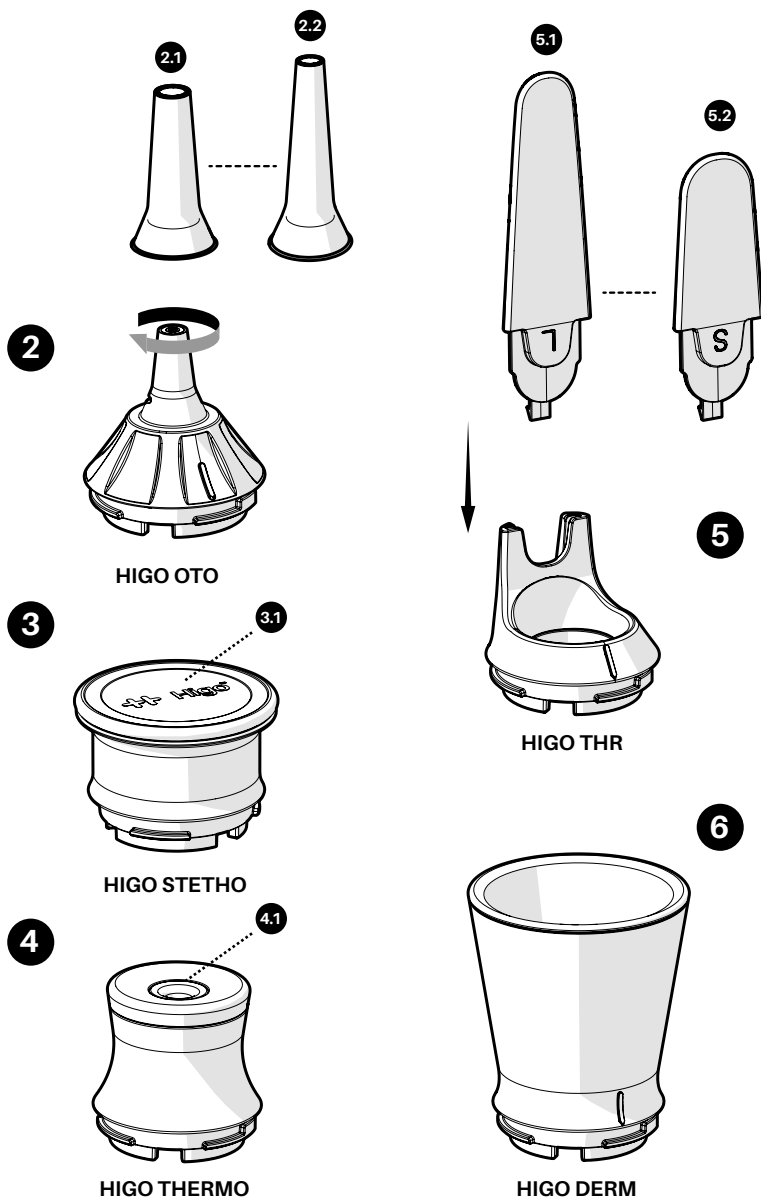
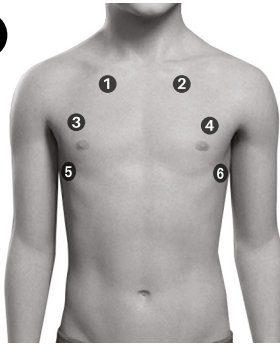
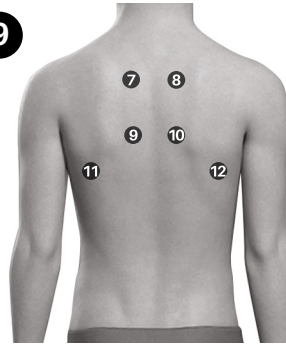
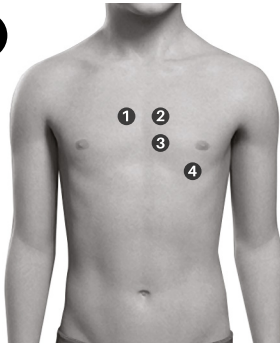
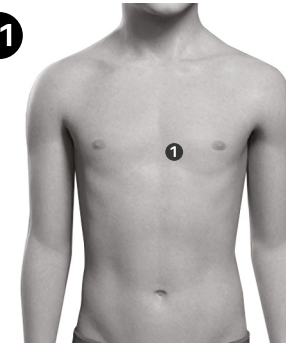
14

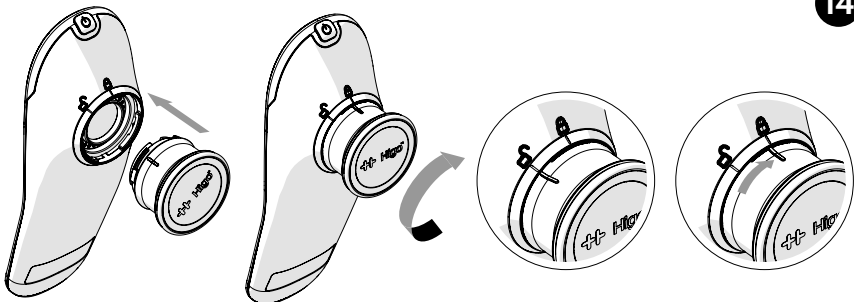
Table of contents

0. Drawings	3
1. Safety information	7
2. Product overview	8
2.1 Definitions	8
2.2 Intended use and indications for use	9
2.2.1 Contraindications	9
2.3 The box components	9
2.4 Overview of the Higo® interface	9
3. How to use the Higo®	11
3.1 Higo® for UNTRAINED OPERATORS - FAMILY OPERATING MODE	11
3.1.1 Prior to first use	11
3.1.2 Using Higo® in FAMILY OPERATING MODE	11
3.2 Higo® for HEALTHCARE PROFESSIONALS - PRO OPERATING MODE	16
3.2.1 Prior to first use	16
3.2.2 Using Higo® in PRO OPERATING MODE	16
3.3 Medical exams	18
3.4 Charging	25
3.5 Cleaning and maintenance	25
3.5.1 Cleaning	25
3.5.2 Calibration of the thermometer	26
3.5.3 Disposable accessories	26
4. Troubleshooting	27
5. Product specification	28
5.1 Specification of the device	28
5.1.1 Wi-Fi network specification	29
5.1.2 Bluetooth Low Energy network specification	29
5.2 Declaration of Conformity	29
5.3 Description of label and packaging marking	32
6. Warranty and service	33
7. Adverse events and medical incidents	33

1. Safety information

Prior to start-up of Higo®, carefully READ THIS USER MANUAL. The proper function and safety operation of this device depends on the user complying with the safety recommendations presented below. Please make sure the device does not appear damaged or broken. If there are any defects please contact the manufacturer.

Warnings: A warning indicates precautions to avoid the possibility of personal injury or death.

Do not use the device if the device packaging, the device itself or any of its accessories are damaged or not intact.

Do not immerse the device or any of its accessories in water or any other liquid. Do not use a wet device. If the Higo® is wet, wipe off all moisture and wait at least 48 hours before operating or contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

The set contains no user-serviceable components. Do not disassemble, service or repair this product or any of its parts. Do not remove the covers of the device.

Before each use the outer surface of the device components should be checked. Do not use parts with rough surfaces, sharp edges or protrusions, which may cause harm.

The eyes should not be exposed to the illumination light of the product during operation.

If the device does not work properly, discontinue its use and contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

Always use only accessories and disposables delivered by HigoSense. Do not use accessories which are not supplied or recommended by the manufacturer. Connect only authorized accessories to the device.

The device includes small parts which may be swallowed. Keep away from children.

Follow the maintenance instructions specifically described in this User Manual.

Do not drop or insert any objects into any device opening.

Do not use the device during defibrillator use. It's not resistant to electrical shock.

Use the Higo® only on intact skin.

Do not push the otoscope with the specula too deep into the ear. Be careful not to touch the eardrum. If the ear canal is blocked by earwax, stop the exam and try again after cleaning the ear canal.

Do not use the device simultaneously with other electronic equipment, to avoid electromagnetic interference with the operation of the device.

Do not use or stack the device near, on or under other electronic equipment, to avoid electromagnetic interference with the operation of the device.

This User Manual should be read before first use. Keep it for your safety.

The device must only be used for the purpose intended by the manufacturer. Arbitrary use of the product which is inconsistent with the User Manual may result in the loss of guarantee rights and claims in the event of damage.

Follow local governing ordinances and recycling instructions regarding disposal or recycling of the device and its components.

Any incorrectness related to the device, which could result in health deterioration or loss of life of the user or patient, should be reported to the manufacturer (HigoSense) and the competent authority of your country.

Cautions - a caution indicates a condition that may lead to damage to equipment or a lower quality treatment.

Check the condition of the device. If the LCD display is broken, contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com. Before use, make sure that the device is clean. If the LCD display is dirty, clean it with a dry towel.

If the device does not turn on, work properly, or the display freezes and cannot be turned off, contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

If the speaker or light indicator does not work properly, contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

Performing an exam is impossible while charging.

Use only the charger which is included in the set by the manufacturer.

Use only a secure Wi-Fi network.

Remember to choose the correct ear (left or right) before the ear exam. Mixing the data can result in a misleading diagnosis.

Discard all the disposables (ear specula, alcohol wipes) after use.

Before each use or after a change of viewing modes/settings, check the view observed through the otoscope module. It should be a live image with a correct image orientation.

The temperature module (the forehead thermometer) is not intended to be used as a contraceptive method.

Always keep the device and its parts in the storage case.

Notes

Operating conditions: 16–40°C / 10–85% humidity non-condensing / 860–1060 hPa

Storage conditions: -10–45°C / 10–85% humidity non-condensing / 860–1060 hPa

It is recommended to recharge the battery after each use. Storing the device with a drained battery may shorten the battery life.

No module should be placed on the base unit while charging.

2. Product overview

The Higo® is a handheld device consisting of a base unit and exchangeable multi-sensor modules, which allows for the exam of the throat, skin, ears, lungs, heart, abdomen, cough and the measurement of body temperature.

The Higo® device can be used by laypersons and healthcare professionals.

2.1 Definitions

SYSTEM ADMINISTRATOR – any individual employed or authorized by the HEALTHCARE PROVIDER, responsible for administration of the medical systems, including managing patient accounts, accesses etc.

MOBILE APP – HIGO MOBILE APPs or any mobile app of a healthcare provider, that can be compatible with the Higo® device

HIGO MOBILE APP or **HIGO PROFESSIONAL MOBILE APP** – all mobile apps provided by HigoSense that can be downloaded from Google Play or App Store

UNTRAINED OPERATOR (LAYPERSON) – any user of the Higo® device who uses it for the purpose of a medical exam or medical data collection, but does not have any medical background or education

MEDICAL QUALITY IMAGES/SOUNDS – images or audio recordings registered by the Higo® device which present a good enough quality and quantity for a valid medical interpretation to be performed by a physician or medical specialist

HEALTHCARE PROVIDER – a health facility providing you with healthcare services based on Higo®

HEALTHCARE PROFESSIONAL – any individual employed or authorized by the HEALTHCARE PROVIDER, having medical knowledge or experience necessary to perform medical exams in PRO OPERATING MODE

RECEIVING SERVER – a server where the data collected by the Higo® device is sent and stored

FAMILY OPERATING MODE – a mode of operation for UNTRAINED OPERATORS at home with the aim to collect medical data and send them to the doctor for interpretation.

PRO OPERATING MODE – a mode of operation for qualified HEALTHCARE PROFESSIONALS in a clinical and ambulatory environment

2.2 Intended use and indications for use

The Higo® is a multi-sensor device used for diagnostic support which digitally records the following medical data:

- images of the mouth, ear canal and skin,
- audio of the lungs, heart, abdomen, and cough,
- body temperature.

The Higo® is intended for use by UNTRAINED OPERATORS at home and by qualified HEALTHCARE PROFESSIONALS in a clinical and ambulatory environment. The product is not intended for self-diagnosis.

The Higo® is a useful tool for collecting medical data on the first symptoms of colds, flu and other viral and bacterial infections of the respiratory tract, middle ear or skin. A patient must be diagnosed by a physician.

2.2.1 Contraindications

There are no known contraindications to the use of the Higo®.

However, special attention should be paid when selecting accessories (tongue depressing tips, ear specula) and performing throat and ear exams on patients under 2 years of age.

2.3 The box components

The Higo® is an active medical device consisting of a base unit **1** powered by a lithium polymer battery and 5 interchangeable modules:

- Otoscope module labelled as **HIGO OTO** **2**
- Stethoscope module labelled as **HIGO STETHO** **3**
- Thermometer module labelled as **HIGO THERMO** **4**
- Throat module labelled as **HIGO THR** **5**
- Dermhood labelled as **HIGO DERM** **6**

Besides the device, the Higo® set includes:

A. Accessories:

- 2 x large exchangeable tongue depressing tips for older children, teenagers and adults **5.1**
- 2 x small exchangeable tongue depressing tips for infants and toddlers **5.2**
- 20 x 4.0 mm disposable ear specula for older children, teenagers and adults **2.1**
- 10 x 2.5 mm disposable ear specula for infants and toddlers **2.2**
- 10 x alcohol wipes.

B. Storage case

C. Medical grade USB-C wall charger with USB Type-C cable Globtek GTM86100-10VV-W2E **7**

D. User Manual

E. Quick Guide

2.4 Overview of the Higo® interface

Note

All drawings can be found in section 0 of this User Manual.

Power button **1.5** – allows you to turn the device on and off.

USB-C charging slot **1.8** – allows you to connect the battery charger.

Audio jack (3.5 mm) 1.9 – allows you to connect headphones. The headphones are used to listen to the sounds generated by the stethoscope module in real time (for detailed requirements regarding headphones, see the Section 5.1).

LED indicator 1.1 – used to inform you about the status of the battery (details are listed below).

LED indicator status	Description
No LED light	The device is ready to be used
Blinking ORANGE	The battery status is below 10%, connect the device to the wall charger
Steady ORANGE	Charging in progress
Steady GREEN light	The device is fully charged
Blinking GREEN light	The device is in sleep mode
Steady RED light	Safety failure mode, see the <i>Troubleshooting Section</i>

Main screen 1.2 :

-  **Wi-Fi icon** – informs about the Wi-Fi connection status and level; if not connected then 
-  **Battery level icon** – informs about the current battery level
-  **Settings menu icon** – clickable icon with the list of device settings menu:
 -  **Wi-Fi** – enables Wi-Fi networks management (adding or removing networks or connecting network from the list)
 -  **Account** – enables user account management (pairing and deactivating the connection with user account)
 -  **Software** – enables verification of firmware version
 -  **Language** – enables changing the language of user interface
 -  **Factory reset** – enables restoring the device to factory settings
-  or **START** button or **REC** button – used to move forward, start an exam or record images/sounds during an exam
-  or **BACK** – used to move backwards
-  or **CANCEL** – used to cancel an exam
- SKIP** button – used to e.g. skip single auscultation points during a heart, lungs or abdomen exam
-  or **REDO** button – used to repeat any exam

Note

The icons presented in this manual are used for explanatory purposes. The actual icons on the device may differ from the images in this manual.

Firmware update

The device firmware will be subject to updates to the newest version. Notifications about the new release of firmware updates will be published in the MOBILE APP and on your Higo® device.

Note

As long as an update has not been completed on your device, you will not be able to perform any exams. For a successful firmware update a Wi-Fi connection is required as well as a device battery level of a minimum of 50%.

Battery saving modes

After detecting the idle state (after 60 seconds), the device turns off the power-hungry process: the LCD display and the camera (with LEDs, if they are turned on). The sleep mode is signalled by the blinking of a green LED ^{1.1}.

If, within the next 10 minutes after entering the mode of limited power consumption,

- the device detects motion, it turns on the display and, if necessary, it restores the camera view,
- the device does not detect any motion, it turns itself off automatically.

Note

In case the device turns off automatically, after 2 hours, your open exam might be automatically closed. It's caused by the fact that medical data becomes obsolete as your health condition may change suddenly.

3. How to use the Higo®

The Higo® device is intended for use by UNTRAINED OPERATORS at home and by qualified HEALTHCARE PROFESSIONALS in a clinical and ambulatory environment. The OPERATING MODES for those two groups of users may slightly differ, which has been explained in the section 3.1 and 3.2.

3.1 Higo® for UNTRAINED OPERATORS - FAMILY OPERATING MODE

3.1.1 Prior to first use

- a. After unpacking the product, check it for completeness. The box should contain all of the items described in Section 2.3.
If you find any of these items is missing, contact your HEALTHCARE PROVIDER or HigoSense immediately - via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.
- b. Charge the device for at least 1 hour before the first use.
- c. Download the MOBILE APP onto your smartphone or tablet from Google Play or App Store.

Note

*A Wi-Fi router connected to the Internet is required to set up the mobile application. **Your Wi-Fi router must be 802.11. b/g/n on a 2.4 GHz frequency.** 5 GHz networks are not supported. For more details about Wi-Fi router compatibility, refer to Section 5.1.1.*

- d. Log in or register to the MOBILE APP using the credentials provided by HigoSense or your HEALTHCARE PROVIDER.
- e. Create profiles of patients. If you use the HIGO MOBILE APP, tap the Add Profile icon on the main screen. Follow the instructions and fill in the required data.

Note

Instructions for creating profiles on the MOBILE APPS, if different from those presented in this User Manual, will be provided by your HEALTHCARE PROVIDER. If you have not received any instructions, please contact your HEALTHCARE PROVIDER or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

- f. Pair the Higo® device with the MOBILE APP (see Section 3.1.2).

3.1.2 Using Higo® in FAMILY OPERATING MODE

Connecting the Higo® device with a Wi-Fi network and pairing with the MOBILE APP

The pairing process between the Higo® and the MOBILE APP happens only once. Please make sure you have the MOBILE APP installed on your mobile. It's required for administrative purposes, e.g. connecting to Wi-Fi or managing patient profiles.

Note

After the factory reset, it is required to repeat the pairing process.

In order to use the Higo® in any new location, you must connect the device with the local Wi-Fi network.

Note

If you change the Wi-Fi password, router or perform any changes to the wireless network, you will need to repeat the process of adding a new Wi-Fi network.

The process of pairing a Higo® device with the HIGO MOBILE APP has been described below.

Note

Instructions for pairing the Higo® device with other MOBILE APPS will be provided by your healthcare provider.

Note

Before pairing, make sure that your mobile phone is connected to a Wi-Fi network and has Bluetooth turned on.

A. Pairing process

Note

Please remember, during the pairing process the Higo® device and the mobile phone have to be close to each other!

- Turn on the Higo® device by pressing the ON/OFF button **1.5** for 1 second.
 - If you have not downloaded the MOBILE APP before, please do it now. You can download it from Google Play or App Store or use a QR code on the device screen.
 - Tap  to move forward.
- Go to the HIGO MOBILE APP and follow the steps below:
 - Log in (or register) to your HIGO MOBILE APP
 - In the app, the connection process has been initiated. If a button to confirm the process appears on the screen, tap to move forward.
 - Choose the Wi-Fi network you want to connect from the list, enter password and click **Connect**.
- The Wi-Fi network will be added and the Higo® will be paired with the HIGO MOBILE APP.
- Data synchronisation will start thanks to active Wi-Fi connection and may take a few minutes. Especially in case of first usage, it might be required to update the device firmware, as described in section 2.4.

B. Connecting to a new (additional) Wi-Fi network

- Turn the Higo® device on by pressing the ON/OFF button **1.5** for 1 second.
 - Tap the settings menu 
 - Choose **Wi-Fi** and tap **Connect new**
- Log in to the HIGO MOBILE APP, follow the steps below:
 - Tap the settings menu 
 - Choose **Higo Device** and tap **Add Wi-Fi**
 - Choose Wi-Fi from the list, enter password and click **Connect**.
- A newly added Wi-Fi network will be added to the device.

For both procedures, data synchronizing will start automatically and take a few minutes. The device has been successfully paired when the Wi-Fi icon changes (on the Higo® device) from  to .

As soon as the Higo® device is successfully connected to Wi-Fi and the MOBILE APP, the list of profiles will be updated to the Higo® device. Then you can start the exam procedure.

Note

The ongoing Wi-Fi connection is required to perform any exam and send the data to the RECEIVING SERVER.

Types of medical exam procedures in the FAMILY OPERATING MODE

There are 2 general types of exams available with the Higo®:

- Quick Check** – singular and individual procedure of registering medical data separately and independently from each other, using a specific module, e.g. temperature.
- Full Medical Exam** – comprehensive set of exams of registering medical data which is performed with the intention of remote interpretation and diagnosis by a physician. Full Medical Exam consists of specific medical exam steps connected to symptoms reported by a patient. The configuration of the exams which should be performed and registered with respect to the symptoms are presented in the Table 1 below (for CHILDREN below 14 years old and TEENAGERS and ADULTS above 14 years old).

Starting an exam

After setting up your device and signing in to the MOBILE APP, everything is ready to start the exam.

- Place the Higo® device in front of you. If the device is off, turn it on by pressing the ON / OFF button  for 1 second.
- Synchronisation with the server will be initiated after a few seconds. If it has not been started automatically after a few minutes, please make sure you have access to the Wi-Fi network you added when pairing the device with the app and the connection is active.
- Make sure the profile of the desired patient is created. If yes, select and tap the profile of the patient you want to examine. If no, log in to the MOBILE APP and see Section 3.1.1.e.

Table 1. Symptoms and exams: CHILDREN below 14 years old and TEENAGERS and ADULTS above 14 years old.

Reportable symptoms	ADULT Set of exams on Higo® device	CHILD Set of exams on Higo® device
Fever (above 38°C)	Temperature measurement Heart Rate Heart auscultation Lungs auscultation	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Throat exam Ear exam Abdomen auscultation
Cough	Temperature measurement Heart Rate Lungs auscultation Heart auscultation Throat exam Cough registration	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Throat exam Ear exam Cough registration Abdomen auscultation

Reportable symptoms	ADULT Set of exams on Higo® device	CHILD Set of exams on Higo® device
Runny nose	Temperature measurement Heart Rate Throat exam	Temperature measurement Heart Rate Lungs auscultation Throat exam Ear exam
Sore throat	Temperature measurement Heart Rate Throat exam	Temperature measurement Heart Rate Lungs auscultation Throat exam Ear exam
Hoarseness	Temperature measurement Heart Rate Throat exam	Temperature measurement Heart Rate Throat exam
Headache	Temperature measurement Heart Rate Heart auscultation	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Throat exam Ear exam
Earache	Temperature measurement Heart Rate Throat exam Ear exam	Temperature measurement Heart Rate Lungs auscultation Throat exam Ear exam
Skin lesions	Temperature measurement Heart Rate Skin exam Lungs auscultation	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Throat exam Skin exam
Swelling (Fluid retentions)	Temperature measurement Heart Rate Lungs auscultation Heart auscultation Skin exam	Temperature measurement Heart Rate Lungs auscultation Heart auscultation Skin exam
Behavior change	n/a	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Ear exam
Shortness of breath	Temperature measurement Heart Rate Lungs auscultation Heart auscultation Throat exam	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Throat exam

Reportable symptoms	ADULT Set of exams on Higo® device	CHILD Set of exams on Higo® device
Diarrhea	Temperature measurement Heart Rate Abdomen auscultation	Temperature measurement Heart Rate Abdomen auscultation Lungs auscultation Throat exam Ear exam
Vomiting / Nausea	Temperature measurement Heart Rate Abdomen auscultation Throat exam	Temperature measurement Heart Rate Abdomen auscultation Lungs auscultation Throat exam Ear exam
Stomach ache and / or testicular pain	Temperature measurement Heart Rate Abdomen auscultation	Temperature measurement Heart Rate Abdomen auscultation Lungs auscultation Throat exam Ear exam
Difficulty / pain when urinating, or a change in urine consistency	Temperature measurement Heart Rate Abdomen auscultation	Temperature measurement Heart Rate Abdomen auscultation Lungs auscultation
Injury or wound	Temperature measurement Heart Rate Skin exam	Temperature measurement Heart Rate Lungs auscultation Skin exam
Loss of smell and/or taste	Temperature measurement Heart Rate Throat exam Lungs auscultation	Temperature measurement Heart Rate Throat exam Ear exam
Muscle and/or joint pain	Temperature measurement Heart Rate Throat exam Heart auscultation Lungs auscultation Contextual picture (of the painful area)	Temperature measurement Heart Rate Throat exam Ear exam Heart auscultation Lungs auscultation Contextual picture (of the painful area)
Excessive fatigue	Temperature measurement Heart Rate Throat exam Heart auscultation Lungs auscultation Contextual movie	n/a
Bug bite / sting	Temperature measurement Heart Rate Lungs auscultation Skin exam	Temperature measurement Heart Rate Lungs auscultation Skin exam
Other	Temperature measurement Heart Rate Lungs auscultation Heart auscultation	Temperature measurement Heart Rate Heart auscultation Lungs auscultation Abdomen auscultation Ear exam

3.2 Higo® for HEALTHCARE PROFESSIONALS - PRO OPERATING MODE

PRO OPERATING MODE is foreseen to be used as a part of integrated healthcare IT system of a HEALTHCARE PROVIDER. For an effective usage of Higo® in PRO OPERATING MODE, a close cooperation with your local SYSTEM ADMINISTRATOR is required, who shall ensure a systematic operator accounts and patients profiles management and update, including delivering login credentials, accounts configuration, pairing and unpairing with Higo® device. For the HEALTHCARE PROFESSIONALS, it's advised to use Higo® in PRO OPERATING MODE in close cooperation with the SYSTEM ADMINISTRATOR designated by the HEALTHCARE PROVIDER or with HigoSense directly.

3.2.1 Prior to first use

- a. After unpacking the product, check it for completeness. The box should contain all of the items described in Section 2.3.
If you find any of these items is missing, contact your SYSTEM ADMINISTRATOR or HigoSense immediately – via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.
- b. Charge the device for at least 1 hour before the first use.
- c. Download the MOBILE APP (e.g. HIGO PROFESSIONAL MOBILE APP) onto your smartphone or tablet from Google Play or App Store.

Note

*A Wi-Fi router connected to the Internet is required to set up the mobile application. **Your Wi-Fi router must be 802.11. b/g/n on a 2.4 GHz frequency.** 5 GHz networks are not supported. For more details about Wi-Fi router compatibility, refer to Section 5.1.1.*

- d. The Higo® device should be paired with the MOBILE APP by your SYSTEM ADMINISTRATOR. In case it is not, please contact your SYSTEM ADMINISTRATOR..

3.2.2 Using Higo® in PRO OPERATING MODE

Connecting the device with operator account and patients' profiles management

The usage of Higo® by HEALTHCARE PROFESSIONALS will be supported by SYSTEM ADMINISTRATOR, who will especially take care of the process of pairing and uploading and updating the list of patients described in section 3.2.1.

- a. Connect the device with your personal operator account following below instructions:
 1. Turn on the Higo® device by pressing the ON/OFF button **1.5** for 1 second.
 - QR code will be displayed on the main screen of the device
 2. Go to the HIGO PROFESSIONAL MOBILE APP:
 - Log in to your operator account in the HIGO PROFESSIONAL MOBILE APP
 - Tap **Connect the device** button
 - Scan the QR code visible on the device
 3. The device is connected with the HIGO PROFESSIONAL MOBILE APP.
 4. Go to the Higo® Device and tap **Synchronize** to initiate the data synchronization. As soon as the Higo® device is successfully connected to Wi-Fi, the list of profiles will be updated into the Higo® device. Then you can start the exam procedure.

Note

Instructions for pairing with MOBILE APP, if different from those presented in this User Manual, will be provided by your SYSTEM ADMINISTRATOR. If you have not received any instructions, please contact your SYSTEM ADMINISTRATOR or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

- b. Patients' profiles are being uploaded during data synchronization on the Higo® device. They are managed by the SYSTEM ADMINISTRATOR.
- c. In case the patient profile you are going to examine does not exist, it should be created manually in the HIGO PROFESSIONAL MOBILE APP. Go to the HIGO PROFESSIONAL MOBILE APP, tap the **Add Patient** icon on the main screen and follow the instructions. Then go to the device and tap **Synchronize**. Please note, in order to perform a synchronisation process, a Wi-Fi connection is required on both device and mobile phone.

Note

Instructions for creating profiles, if different from those presented in this User Manual, will be provided by your SYSTEM ADMINISTRATOR. If you have not received any instructions, please contact your SYSTEM ADMINISTRATOR or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

Types of medical exam procedures in the PRO OPERATING MODE

The PRO OPERATING MODE allows one type of exam procedure:

1. **Medical Exam Procedure** – comprehensive set of exams of registering medical data which is performed with the intention of medical interpretation and diagnosis by a physician. Medical Exam Procedure consists of exams determined and decided by a HEALTHCARE PROFESSIONAL.

Performing Medical Exam Procedure

- a. Select a patient profile and confirm your choice
- b. Start the **Medical Exam**
- c. Perform necessary exams from the list by tapping a particular exam and follow the instructions on the screen.

Note

In the PRO OPERATING MODE it is an operator (HEALTHCARE PROFESSIONAL) who decides about a set of exams to be performed during the Medical Exam Procedure based on their medical expertise. It is recommended, however, to follow the combination of medical exams for reported symptoms as presented in the table in section 3.1.2.

- d. As soon as all required exams are completed, tap **Finish the exam**.
- e. The registered data are stored in the device's internal memory until next synchronisation.
- f. In order to synchronize stored data: go to the device and tap **Synchronize**. Please note, in order to perform a synchronisation process, a Wi-Fi connection is required on both device and mobile phone.

Note

For successful synchronization it is required to have established a stable Wi-Fi connection. In case the Wi-Fi connection is not available the registered exams will be saved in the internal memory of the device. Remember to synchronize completed exams as soon as the Wi-Fi connection is available.

3.3 Medical exams

Before each use, disinfect and clean the surface of the Higo® device, modules and tongue depressing tips using 70% isopropyl alcohol wipes or cloth/cotton swabs dipped in rubbing alcohol. You can use the alcohol pads attached to the device set. Make sure that the lenses located on the front of the otoscope module and the camera on the base unit are clean.

If any symptoms present are severe and/or long-lasting please contact a primary care physician or emergency medical services. The Higo® is not designed for self or automatic diagnosis and to be used in life-threatening situations. Failure to contact appropriate medical authorities may result in serious bodily injury or death.

In case that a patient's skin surface is injured, burned, cut, soaking, with abrasions, with open sores or represents any other wounds in a place where the exam expects a physical contact with a device or any accessories, the exam should be withdrawn.

Note

Please remember, the device does not allow the review of any data. The data is temporarily stored on the internal SD CARD before it will be sent to the RECEIVING SERVER.

Note

Always verify the quality of the Wi-Fi connection before starting the exam. If the Wi-Fi indicator shows poor quality (below 2 bars), there may be a problem with data transfer.

A. Temperature measurement

The purpose of the temperature exam is to measure the body core temperature by scanning the forehead and especially the temporal artery which is located on the left and right temple.

The measurement takes place by registering temperature values with a high frequency by an advanced infrared sensor with the aim of scanning the temporal artery which remains a reference for the body core temperature and uploads it to the RECEIVING SERVER, thanks to which they can be later interpreted by a physician or stored for medical comparison using a medical data viewer.

Before taking the TEMPERATURE MEASUREMENT, make sure that:

- The patient and the device remain in the same ambient temperature for at least 10 minutes,
- The skin of the forehead is dry,
- The measurement is not performed within at least 30 minutes of showering, bathing, exercising or any other physical activity which could influence the forehead surface temperature,
- You do not use the thermometer on scar tissue, open sores or abrasions,

- The thermometer is kept in a storage case in ambient temperature. If the thermometer module was stored in a sunny or very cool place (e.g. next to a radiator, on a windowsill, in the boot of a car), it should be allowed to acclimatize to room temperature for 20 minutes before the measurement is carried out.
- You do not take multiple scans in rapid succession, as these will change the temperature of the skin and results variability of the readings. It is important to wait 10 seconds before repeating the scan to allow the skin to recover.

Step 1: Select the temperature measurement from the exams list.

Step 2: Connect the temperature module **4** to the base unit **1** as presented at Fig. **14**.

Step 3: Read the instructions and information presented on the device screen carefully and click **>** to take a measurement.

Step 4: Place the thermometer module on the centre of the forehead above the line of the eyebrows and click START to begin the measurement. Move the thermometer along the forehead towards the ear tip. Make sure that the temporal artery on the temple is covered during this movement.

The measurement takes 5 seconds and stops automatically with a soft sound. The measurement result is presented on the screen in Celsius degrees.

You can measure the left or right side of the forehead.

In order to obtain higher accuracy, repeat the measurements at least 2 times on the left and 2 times on the right side of the forehead.

Due to biological variety, results from the left and right side will most probably vary from each other by up to 1°C which is completely normal. **Remember to always consider the highest result.**

Age	Range of normal temperature	Low-grade fever	High fever
0-3 months	36-38°C	-	>38°C
3-36 months	36-38°C	38-39°C	>39°C
36+ months	36-38°C	38-39.5°C	>39.5°C

Step 5: Click **>** to save and send the data to the RECEIVING SERVER or REDO to start the measurement once again. You can do as many measurements as you want.

You can do the measurement either with the physical contact of the thermometer with the skin or contactless at 1 cm above the skin. Remember to always keep the same distance from the skin along the contactless measurement in order to ensure reliable results.

Keep the thermometer module clean and protect it from dust or other small contaminants for high accuracy and repeatability of measurements. In the event of contamination inside the aluminium cone or the temperature sensor **4.1**, use a water moistened cotton swab and clean it gently. Do not use body or dish soap or detergents.

Note

The results may vary between measurements even on the same person on the right and left side of the forehead. It depends on multiple factors such as ambient temperature and humidity, medication taken, physical condition, menstrual cycle, and other factors.

B. Auscultation

Heart Rate, Heart exam, Lungs exam, Abdominal exam

There are 4 auscultation exams available to perform using the Higo® device:

1. Heart rate
2. Heart auscultation exam
3. Lungs auscultation exam
4. Abdominal auscultation exam

The purpose of the auscultation exams is to record the sounds from auscultation points on the chest, back or abdomen in such a way that a physician is able to listen and interpret them via headphones.

The device registers audio sounds using the stethoscope module equipped with a membrane similar to a classic stethoscope. The audio recordings are uploaded to the RECEIVING SERVER thanks to which they can be later interpreted by a qualified medical professional or stored for medical comparison using a medical data viewer.

For all auscultation exams it is possible to listen to the sounds while examining using headphones plugged in via a mini-jack slot placed in the bottom of the device.

Before starting the AUSCULTATION EXAM, make sure that:

- There is no ambient noise, like music, speaking or other noticeable sounds in the background which could affect the auscultation sounds,
- The patient is undressed from the waist up. **No bra is allowed during the auscultation exams (lungs and heart),**
- The stethoscope membrane is fully applied to the body with its whole surface,
- You are not placing the stethoscope module on the sternum,
- You are not scratching, tapping or scraping the enclosure of the device which can generate internal noises.

Note

The heart rate readout represents an estimation of the physiological heart rate. The readout is based on the same acoustic signal that is delivered to the physician for the diagnostics purposes.

Note

The stethoscope module does not provide diagnostic information, but it is a diagnostic aid. It amplifies body sounds, which the physician must interpret to make a diagnosis.

Note

Stethoscope does not provide vital signs monitoring and alarms for potential life threatening situations in which medical intervention is necessary.

Step 1: Select the Heart Rate, Heart, Lungs or Abdominal exam from the menu on the device screen.

Step 2: Connect the stethoscope module **3** to the main unit **1** as presented at Fig. **14**.

Step 3: Read the instructions and information presented on the screen carefully and click **>** to start the exam.

Step 4: Instruct the patient to sit straight and start taking deep breaths whilst doing the Lungs exam: breathe in and out through the mouth. For the Heart Rate, Heart and Abdominal exam instruct the patient to breathe normally.

Step 5: Place the stethoscope module on the naked skin firmly on the auscultation points presented on the screen and in the figures in section 0, in the correct order, ensuring full contact of the stethoscope membrane (Fig. 3.1) with the skin surface.

Heart Rate auscultation point is presented on the picture 11 in section 0.

Heart auscultation points are presented on the picture 10 in section 0.

Lungs auscultation points are presented on the pictures 8 and 9 in section 0.

Abdominal auscultation points are presented on the picture 12 in section 0.

The exam and data registering starts automatically at each auscultation point only when contact between the stethoscope membrane and skin is established.

The sound registering should stop automatically after ca. 10–60 seconds as soon as the MEDICAL QUALITY SOUND is recognized by the device.

In case the MEDICAL QUALITY SOUND is not recognized within 60 seconds, the recording stops automatically and information about the potential reason will be shown on the screen. There may be 2 general reasons for not recognizing the MEDICAL QUALITY SOUNDS:

1. Low quality of the sounds and/or MEDICAL QUALITY SOUND is hardly audible
 - a. During the Lungs exam, the patient is not taking deep breaths
 - b. The stethoscope module was applied in an incorrect position on the body
2. Too much ambient noise in the background

In such cases you are advised to REDO the registration, making sure that you're following the instructions carefully, or ACCEPT the registration anyway to save and send the data to the RECEIVING SERVER without medical quality confirmation.

In the event of losing contact between the stethoscope and skin surface before the exam is stopped automatically, it is necessary to REDO the recording.

Note

*It is acceptable for small children to lie down on their back or front during the auscultation exams in order to stay calm and cooperate with the operator.
During the Lungs exam, the time required for automatic registering stops may vary between patients. It may depend on age, physical condition, fat tissue and medical history.*

Cough registration

Cough registration is a supporting exam which helps in the interpretation of a cough type and origin.

The purpose of the cough registration is to record the cough sound using the stethoscope module in such a way that a physician is able to listen to and interpret it via headphones.

Before starting the COUGH REGISTRATION, make sure that:

- There is no ambient noise, like music, speaking or other noticeable sounds in the background which could affect the auscultation sounds,
- You are not scratching, tapping or scraping the enclosure of the device which can generate internal noises.

Step 1: Select the Cough exam from the menu on the device screen.

Step 2: Connect the stethoscope module 3 to the main unit 1 as presented at Fig. 14.

Step 3: Read the instructions and information presented on the screen carefully and click > to start the exam.

Step 4: Place the stethoscope module ca. 5 cm from the patient's mouth (not directly on the stethoscope membrane) as presented at Fig. 13 and click the START button on the screen.

Step 5: Instruct the patient to cough loudly 2–3 times, as he/she would cough into a microphone.

Step 6: Click STOP to terminate the recording, otherwise the registering stops automatically after 10 seconds.

Step 7: Click ➤ to accept the recording and save and send the data to the RECEIVING SERVER or click REDO to repeat the exam.

C. Left and right ear exam

The purpose of the ear exam is to depict the tympanic membrane in such a way that a physician is able to interpret its physical condition and healthiness by judging anatomy, size, colour and surface. The device registers images with a high frequency using an embedded camera and the otoscope module and uploads them to the RECEIVING SERVER thanks to which they can be later interpreted by a qualified medical professional or stored for medical comparison using a medical data viewer.

Before starting the EAR EXAM, make sure that:

- You use a fresh and clean disposable ear speculum,
- You have selected the correct size of speculum with respect to the age of patient,
- You're able to recognize the left and right ear of the patient,
- You do not use a disposable speculum more than once. Throw it away after the exam of the patient,
- You DO NOT attempt the exam without a speculum applied,
- You DO NOT look directly into the otoscope tip when the light is on.

Step 1: Connect the otoscope module ② to the main unit ① as presented at Fig. 14.

Step 2: Apply the correct speculum size according to the patient's age, by inserting the speculum onto the otoscope tip and twisting it until it is locked into place.

Choose the large speculum marked with size "4.0 mm" 2.1 for patients above 3 years of age.

Choose the small speculum marked with size "2.5 mm" 2.2 for patients below 3 years of age.

Step 3: Read the instructions and information presented on the screen carefully and click ➤ to start the exam.

Step 4: Grab the device sideways in your left hand whilst examining the left ear or in the right hand whilst examining the right ear and slowly insert the otoscope tip with the speculum into the ear canal.

Important

Before inserting the otoscope into the ear canal, grab the earlobe tightly and pull firmly to the back or diagonally in order to straighten the ear canal. Thanks to this, the ear exam will be more comfortable for the patient and more efficient for the operator.

Step 5: Search for the optimal view of the tympanic membrane by moving the otoscope tip gently forwards and sideways inside the ear canal.

Registering starts automatically as soon as the camera review is activated on the screen. The device recognizes the MEDICAL QUALITY IMAGES of the tympanic membrane and confirms it showing a green symbol.

Step 6: Click the "STOP" button on the screen to complete the exam or wait 30 seconds so the exam is finished automatically. An image of the best MEDICAL QUALITY will be displayed on the screen, or information about no proper image recorded.

Click ➤ to accept the exam and send it to the RECEIVING SERVER.

In case there is no MEDICAL QUALITY IMAGE recorded, click REPEAT to redo the exam or "SEND CURRENT EXAM" in order to send it to the RECEIVING SERVER without confirmed MEDICAL QUALITY.

The MEDICAL QUALITY IMAGE of a tympanic membrane contains clearly visible:

- umbo
- manubrium of malleus
- light reflex (on the healthy tympanic membrane)

Note

Due to naturally occurring variability in anatomical structures of ear canal and tympanic membranes, and the effects of illness or past surgical interventions, it may happen that not all structural elements as well as the light reflex are visible in all examined subjects. However, the overlying principle of a tympanic membrane image exam is that the collected data will portray a current and accurate image of an individual subject's tympanic membrane.

Pay special attention to earwax, which is the most common reason preventing the registering of MEDICAL QUALITY data.

Remember to start registering before inserting the otoscope tip into the ear canal to avoid clicking the screen during the exam, which can cause discomfort to the patient and extend the exam unnecessarily.

Note

Earwax is a natural product of your body. Its presence in the ear canal is normal until it causes hearing problems or pain. If there is too much earwax in the ear canal, making it impossible to register for the medical quality exam, try to clean the ear canal or visit an ENT doctor.

D. Throat exam

The purpose of the throat exam is to depict the oral cavity in such a way that a physician is able to interpret its physical condition and healthiness by judging anatomy, size, colour and surface. The device registers consecutive images using an embedded camera and illumination system and uploads them to the RECEIVING SERVER thanks to which they can be later interpreted by a physician or stored for medical comparison using a medical data viewer.

Before starting the THROAT EXAM, make sure that:

- The tongue depressing tip is clean and disinfected properly,
- It is at least 15 minutes after the patient's last meal,
- You've selected the correct tongue depressing tip size with respect to the age of the patient.

Step 1: Connect the throat module **5** to the main unit **1** as presented at Fig. **14**.

Step 2: Connect the correct tongue depressing tip according to the patient's age, by placing the tip into the throat module slot until it clicks and is locked into place automatically.

Choose the large tongue depressing tip marked with the letter "L" **5.1** for patients above 5 years old.

Choose the small tongue depressing tip marked with the letter "S" **5.2** for patients under 5 years old.

For professional use it is strongly recommended to use wooden disposable tongue depressors placing them into a slot in the throat module the same way as with reusable ones

Step 3: Read the instructions and information presented on the screen carefully and click **>** to start the exam.

Step 4: Place the depressing tip on the tongue to $\frac{3}{4}$ of its length, parallel to the palate and press against the tongue gently and firmly. Ask the patient to say "Ahhh!" without puffing which can cause the camera to fog.

Registering starts automatically as soon as the camera review is activated on the screen. The device recognizes the MEDICAL QUALITY IMAGES of the throat and confirms it showing a green symbol.

Step 5: Click the "STOP" button on the screen to complete the exam or wait 30 seconds so the exam is finished automatically. An image of the best MEDICAL QUALITY will be displayed on the screen, or information about no proper image recorded.

Click  to accept the exam and send it to the RECEIVING SERVER.

In case there is no MEDICAL QUALITY IMAGE recorded, click REPEAT to redo the exam or “SEND CURRENT EXAM” in order to send it to the RECEIVING SERVER without confirmed MEDICAL QUALITY.

The MEDICAL QUALITY IMAGE of a throat contains clearly visible:

- uvula
- tonsils
- back of the throat
- palatine arches

Note

Due to naturally occurring variability in anatomical structures of the oral cavity, and the effects of illness or past surgical interventions (e.g. tonsillectomy), it may happen that not all structural elements are visible in all examined subjects. However, the overlying principle of a throat image exam is that the collected data will portray a current and accurate image of an individual subject's throat.

Remember to start registering before applying the tongue depressing tip into the oral cavity to avoid clicking the screen during the exam which can cause discomfort to the patient and extend the exam unnecessarily.

Note

There are people who cannot stand having their throat examined with the tongue depressing tip. In such cases it is permitted to carry out the exam without the depressing tip attached, however the recognition of MEDICAL QUALITY is less probable.

E. Skin exam

The purpose of the skin exam is to depict skin lesions in such a way that a physician is able to interpret the skin's physical condition. The device registers consecutive images using an embedded camera and illumination system equipped with the dermhood, and uploads them to the RECEIVING SERVER thanks to which they can be later interpreted by a physician or stored for medical comparison using a medical data viewer.

Before starting the SKIN EXAM:

- Make sure that the skin surface is clean, dry, naked. No skin creams, ointments, lotions, emulsions, balsams is allowed to be applied before and during skin exam
- Disinfect the edges of the dermhood before each use

Do not perform the skin exam during a medical treatment with photoactive substances, due to white light emission by the device's illumination system.

Step 1: Connect the dermhood **6** on the main unit **1** as presented at Fig. **14**.

Step 2: Read the instructions and information presented on the screen carefully and click  to start the exam.

Step 3: Apply the dermhood to the skin in such a way that it touches the skin surface. Do not press it so that you make a bulge in the skin.

Step 4: Start registering by clicking START on the screen and slowly move the dermhood over the skin surface in order to depict the whole area of the skin change. If the skin change or lesion is extensive, make sure you have covered the entire area of the change or lesion as well as borders with unchanged healthy skin.

Step 5: Click STOP to terminate the recording, otherwise the registering stops automatically after 30 seconds.

Step 6: Click  to accept the image and send it to the RECEIVING SERVER or click REDO to repeat the exam.

Reviewing recorded exams

On condition that the registered data has been sent to RECEIVING SERVER, it should be possible to review and present data on the MOBILE APP for your own purposes or share it with a healthcare professional for medical assessment and diagnostics.

3.4 Charging

The Higo® has a rechargeable, lithium-ion polymer battery. You can check the battery status by the icon indicator at the top bar on the device's screen.

To recharge the battery, remove the rubber cover **1.4** situated on the bottom of the device and connect the charger **7** cable to the charging outlet **1.8**. The LED indicator will light up steady orange.

For safety reasons, performing an exam is impossible while the device is charging.

Note

It is recommended to recharge the battery after each use. Storing the device with a drained battery may shorten the battery life.

Note

After storing the device for a long period of time without charging, it may not be possible to immediately use your Higo® device due to the discharged battery. In this case, connect your device to the charging. It's recommended to charge the battery regularly in case the device is not used for longer than 6 months.

3.5 Cleaning and maintenance

Before each exam disinfect and clean the surface of the device, modules and tongue depressing tips. By disinfecting you make sure that the surface is free of pathogenic microorganisms such as bacteria, viruses and fungi.

3.5.1 Cleanig

- Turn off the device before cleaning.
- Start with the base unit and modules. First, wipe off any visible contaminants from the device's surface with a damp cloth. Next, allow the device to air dry for at least two minutes. Finally, clean the entire device housing and camera glass cover thoroughly using the provided 70% isopropyl alcohol wipes for at least 30 seconds. After cleaning, allow the device to air dry for at least two minutes in a well ventilated area before performing any exams.
- The most delicate parts of the device to pay a special attention during cleaning are the following:
 - the lenses located on the front of the otoscope module,
 - the cone and sensor of the thermometer module,
 - the camera glass cover.

These parts must be clean and intact to ensure proper results can be obtained. Wipe their surface very gently with a cotton swab moistened with alcohol wipes for at least 30 seconds. Hold the modules upside down to prevent excess moisture from entering the lenses and sensor area. After the alcohol has dried out (minimum after 2 minutes), you can perform an exam.

- d. Clean the reusable tongue depressing tip with soap under running water for at least 30 seconds, and dry well. Next, thoroughly clean the surface of the tongue depressing tip with 70% isopropyl alcohol wipes for at least 30 seconds. Leave to air dry for at least 2 minutes before using the depressing tip for a throat exam.

For cleaning particularly hard-to-reach places, like the module mechanical interface on the base unit, use a cotton swab dipped in minimum 70% isopropyl alcohol for at least 30 seconds. Make sure you allow the device to dry for at least two minutes in a well ventilated area before performing any exams.

Note

The alcohol wipes and ear specula included in the Higo® set are disposable. Always throw away a used tissue or specula immediately after use.

Do not put any objects into the device openings.

Pay a special attention to the module connection mechanism area while cleaning to prevent bacteria and viruses multiplication.

While not in use, the device should be stored in the storage case in a clean and safe location. The device should not be directly exposed to sunlight. Disinfect the surface with alcohol wipes (70% isopropyl). It is recommended that you clean it before each use.

The Higo® device is not intended to be sterilized.

The device and disposables must be kept at a temperature between -10°C and +45°C. Do not expose the device to UV radiation (e.g. from UV lamps used for air disinfection or from the sun).

The device must not be washed or immersed in water.

The set contains no user-serviceable components. Do not disassemble, service or repair this product or any of its parts. Do not remove the covers of the device

3.5.2 Calibration of the thermometer

The thermometer module **4** was calibrated at the time of manufacture. If used according to the instructions, periodic calibration is not required.

The product has a lifetime of 4 years. To maintain functional product performance and comfort of use, after the first 2 years of use, the device should be delivered to the manufacturer's service centre for a technical review, during which the overall condition of the device will be assessed.

3.5.3 Disposable accessories

Note

It's required to use only disposable accessories recommended by the manufacturer and included in the Higo® set. In order to purchase disposable accessories for the Higo® device, please contact your HEALTHCARE PROVIDER or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

4.Troubleshooting

Situation	Solution
The device is not turning on.	Make sure the battery is charged by connecting the device to the charging. Current battery level will be displayed on the LCD screen in the right upper corner. If the battery level is below 10%, the device needs to be charged first. If the battery is full, disconnect the charger and try to turn on the device. If the device still does not turn on, please make sure you are operating under the acceptable conditions described in this manual (temperature, humidity, atmospheric pressure). If the issue persists, please contact your healthcare provider or HigoSense support.
I cannot add the Wi-Fi network.	Make sure the Wi-Fi name and password you are entering match and are correct. Please remember, 5 GHz band is not supported. If the name and password are correct, please remove all the networks from the device, turn it off and on and retry the process of adding Wi-Fi. If the issue persists, please contact your healthcare provider or HigoSense support.
I have added the Wi-Fi network but it's not connecting automatically.	Please try to restart the device. Every time it is turned on, the Wi-Fi connection is initiated again. If the issue persists, please contact your healthcare provider or HigoSense support.
I cannot pair the device with the app.	Make sure that on your smartphone or tablet the Bluetooth is turned on and visible. Please check if the Wi-Fi connection is stable and you are logged in to the correct account. If the issue persists, please turn the device off and on. If it doesn't help, please contact your healthcare provider to ensure your account is active.
The device does not charge.	Please check if you are using the charger provided by the manufacturer. Make sure it's plugged in. Please refrain from using any other charger than the one provided and recommended by the manufacturer. If the issue persists, please contact your healthcare provider or HigoSense support.
The device has been dropped on the floor and does not turn on.	Please contact your healthcare provider or HigoSense support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com .
Although I have connected the module, the device still shows me the instructions to do this.	Make sure that the module you have connected is the right one, i.e. for the chosen exam. If so, disconnect the module, turn the device off and on and retry. If you've chosen the incorrect exam, please close it and choose the correct one. It may also mean that the module has been wrongly connected or there is a technical failure. In both cases, detach the module and turn the device off and on. If the problem persists, please contact your healthcare provider or HigoSense support.
Sending the data to the server takes a very long time.	Check your Wi-Fi connection – if it's too weak, synchronization will take longer. If the synchronization fails, please check in the device settings if the account is correctly paired and the email address is visible there.
The measured temperature is inaccurate or permanently out of range due to small contaminants present on the temperature sensor 4.1.	Using a water moistened cotton swab, clean the sensor gently and make sure that the surface of the sensor is clean and free from small contaminants.
The temperature taken is lower or higher than typical human temperature range (34–43°C).	Please try to redo the exam. If the situation persists, please turn the device off and on. If the issue persists, please contact your healthcare provider or HigoSense support.
The battery led indicator light is steady red.	Please categorically refrain from using the device. Please contact Higo support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com .
The screen is frozen and does not react.	Please try to restart the device. If the problem persists, please contact your healthcare provider or HigoSense support.
If none of the above mentioned actions solved the problem, please try to perform Factory Reset. It brings the device to the basic state which may solve some issues. Please note, after Factory Reset, you have to pair the device with your account and add Wi-Fi network.	

5. Product specification

5.1 Specification of the device

Dimensions	180 mm x 61 mm x 52 mm
Weight	165 g
Display	2.8"
Audio output port	3.5 mm standard headphone connector
Power supply	USB-C plug, wall charger, USB-C cable input: AC 100-240V, 50/60Hz, 0.3A output: DC 5V, 2A
Wireless connection frequency band	2.4 GHz (see Section 5.1.1)
Base unit	
Max. image resolution	1920x1080 px
Illumination	8 white SMT LEDs ring (5700K)
Stethoscope module	
Frequency range	5–5000 Hz
Thermometer module	
Operating temperature range	34–43°C
Laboratory accuracy	+/- 0.3°C
Operating conditions	
Temperature	16–40°C
Humidity	10–85% non-condensing
Atmospheric pressure	860–1060 hPa (battery)
Storage conditions	
Temperature	-10–45°C
Humidity	10–85% non-condensing
Atmospheric pressure	860–1060 hPa (battery)
Battery life	2 years
Warranty	2 years
Additional accessories included	Storage case USB-C wall charger 20 x 4.0 mm specula 10 x 2.5 mm specula 2 x large tongue depressing tips 2 x small tongue depressing tips 10 x alcohol wipes
Requirements for device to install MOBILE APP	Smartphone or tablet with a Wi-Fi and Bluetooth module - RAM min. 100 MB
Requirements for headphones	Over-ear headphones 3.5 mm mini-jack stereo Cable with min. 1 m Impedance: 32 Ohm Frequency response: 20 Hz – 20 kHz

5.1.1 Wi-Fi network specification

The following Wi-Fi network specifications are supported by the Higo® device:

1. IEEE 802.11n standard is recommended for optimal performance.

IEEE 802.11g standard is also supported.

Radio operates in 2.4 GHz band only.

Please note that 5GHz / AC standard is not supported.

2. Transmitter technical characteristics:

Frequency band: 2400 MHz – 2483.5 MHz.

Maximum output power: 20 dBm (100 mW).

3. Antenna information:

Type: Internal 2.4 GHz.

Gain: 0.712 dBi.

Wi-Fi access points have to support the above technologies and work in the same frequencies range. Failure to provide the above network characteristics may lead to an inability to establish a network connection or Wi-Fi connection breaks, thus limiting the possibility of the device to connect and transmit exam data.

The Higo® device should be connected to a Wi-Fi network that is made exclusively for it, to minimize the potential risks of other equipment's influence. This influence may present additional risks of unknown level for the operator, patients or third parties. It is the responsibility of the user to identify, analyse, evaluate and control those risks. Any change introduced in the network configuration may lead to new risks and should be evaluated to determine its safety. Changes that lead to the re-evaluation requirement include:

- change of Wi-Fi network configuration,
- connecting new devices (clients) to the network, including hardware upgrade of existing devices,
- updating software of the equipment already present in the network. This includes all clients, access points and routers that establish wireless connections.

5.1.2 Bluetooth Low Energy network specification

1. Module information:

Bluetooth radio technology: Bluetooth Low Energy, LE + 2LE

2. Transmitter technical characteristics:

Frequency band: 2400 MHz – 2483.5 MHz.

Maximum output power: 10 dBm (10 mW).

3. Antenna information:

Type: Internal 2.4 GHz.

Gain: 0.712 dBi.

5.2 Declaration of Conformity

The safety, reliability and performance of this device can only be assured under the following conditions:

- a. The device has been used for its indicated use and according to the operating instructions in this User Manual.
- b. All fittings, extensions, readjustments, changes, or repairs have been carried out by HigoSense.

HigoSense declares that the Higo® conforms with the essential requirements and other relevant provisions of the Directive 93/42/EEC and the Directive 2011/65/EU.

The device complies with the requirements of IEC/EN 60601-1 for general requirements for safety of medical electrical equipment and IEC/EN 60601-1-11, IEC/EN 60601-2-18, IEC/EN 62471 and ISO/EN ISO 80601-2-56 for electrical and photobiological safety.

The Higo® has IEC 60601-1-2 Class B compliance which applies to the basic safety and essential performance of medical electrical equipment and ME systems in the presence of electromagnetic disturbances and to electromagnetic disturbances emitted by medical electrical equipment and ME systems.

Declaration - Electromagnetic Emission

The Higo® is intended for use in the electromagnetic environment as indicated below. The user shall ensure that the device is used in such an environment.

Emission test	Compliance	Electromagnetic environment - guidance
RF emission CISPR 11	Group 1	The Higo® generates a radio frequency electromagnetic field for the purpose of operating internal functions only. The emission of radio waves is therefore very low and the likelihood of interference with nearby electronic equipment is low.
RF emission CISPR 11	Class B	The Higo® generates a radio frequency electromagnetic field for the purpose of operating internal functions only. The emission of radio waves is therefore very low and the likelihood of interference with nearby electronic equipment is low.
Harmonic current emissions IEC 61000-3-2	Class A	The Higo® device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Voltage changes, voltage fluctuations and flicker emissions IEC 61000-3-3	Complies	

ESSENTIAL PERFORMANCES

The essential performances of the Higo® are:

- The ability to display live view rather than recorded during ear exam
- The ability to display images in proper orientation
- The ability to measure the body temperature with laboratory accuracy
- The ability to register medically suitable tympanic membrane images

Use of accessories, transducers and cables other than those provided or specified by the manufacturer could result in increased electromagnetic emission or decrease electromagnetic immunity of this equipment and result in improper operation.

Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) from any part of the Higo®, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Declaration of electromagnetic conformity

Immunity test	IEC 60601 level	Compliance level	Electromagnetic environment
Electrostatic discharge (ESD) IEC 61000-4-2	8kV contact 2,4,8,15kV air	8kV contact 2,4,8,15kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for lines power supply	± 2 kV for lines power supply	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5 and 1 kV in differential mode ± 0.5, 1 and 2 kV in common mode	± 0.5 and 1 kV in differential mode ± 0.5, 1 and 2 kV in common mode	Mains power quality should be that of a typical commercial or hospital environment.
Rated power frequency magnetic fields	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply lines IEC 61000-4-11	0% U_T , 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T , 1 cycle; 70% U_T ; 25 cycles at 50Hz and 30 cycles at 60Hz Single phase at 0° Voltage interruptions; 0% U_T ; 250 cycles at 50Hz and 300 cycles at 60Hz	0% U_T , 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T , 1 cycle; 70% U_T ; 25 cycles at 50Hz and 30 cycles at 60Hz Single phase at 0° Voltage interruptions; 0% U_T ; 250 cycles at 50Hz and 300 cycles at 60Hz	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply.
Conducted RF IEC 61000-4-6	3V 0,15 MHz-80 MHz 6V in ISM and amateur radio bands between 0,15MHz and 80 MHz 80% AM at 1kHz	3V 0,15 MHz-80 MHz 6V in ISM and amateur radio bands between 0,15MHz and 80 MHz 80% AM at 1kHz	Portable and mobile RF communications equipment should be used no closer to any part of the Higo®, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF EM fields IEC 61000-4-3	10V/m 80MHz -2,7 GHz 80% AM at 1kHz	3V/m 80MHz -2,7 GHz 80% AM at 1kHz 10V/m 80MHz -2,7 GHz 80% AM at 1kHz	Recommended separation distance: d=1.2150 kHz to 80 MHz d=1.2 80 MHz to 800 MHz d=2.3 800 MHz to 2.5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by electromagnetic site surveys, should be less than the compliance level in each frequency range ⁹ .

Immunity test	IEC 60601 level	Compliance level	Electromagnetic environment
Voltage dips, short interruptions and voltage variations on power supply lines IEC 61000-4-11	9V/m 710MHz, 745MHz, 780MHz, 5240 MHz, 5500 MHz, 5785 MHz	9V/m 710MHz, 745MHz, 780MHz, 5240 MHz, 5500MHz, 5785 MHz	a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device. b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m
	27V/m 385MHz, 28V/m 450MHz, 810 MHz, 870 MHz, 930 MHz, 1720 MHz, 1845MHz, 1970 MHz, 2450 MHz.	27V/m 385MHz, 28V/m 450MHz, 810 MHz, 870 MHz, 930 MHz, 1720 MHz, 1845 MHz, 1970 MHz, 2450 MHz.	

If you would like to receive the declaration of conformity, contact the manufacturer support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com.

Radio frequency and EMC statement

The Higo® device is intended for use in an electromagnetic environment with controlled radio frequency interference. The customer or the user of the device may help prevent electromagnetic interference by maintaining the minimum required distance between the mobile device and cellular radio devices (transmitters). This equipment is not subject to protection from harmful interference and may not cause interference with duly authorized systems. HigoSense is not responsible for any radio or communication interference caused by using an other than recommended charger and battery or by unauthorized modifications of this equipment.

5.3 Description of label and packaging marking

Symbol	Description
 HIGONSENSE SP. Z O.O. ZA JĘCZA 15 00-351 WARSZAWA POLAND	Manufacturer
 2274	CE marking that indicates that a product complies with the Directive 93/42/EEC and the Directive 2011/65/EU. Certified product. Notified body number: 2274
	Serial number identifier (Numbers 4 to 7 indicate the year and week of manufacture of the component. The same numbering is used for the unique identification number of the medical device and the marking plate)
	Lot or batch identifier (numbers from 4 to 7 indicate the year and week of the year of manufacture of the component)
	Model or type reference
IP22	Degrees of protection provided by enclosures (IP Code). The device is protected against solid objects over 12.5 mm (e.g. a finger) and protected against falling drops of water, if the case is disposed of at any angle up to 15 degrees from vertical.

Symbol	Description
	Direct current
	Type BF applied part
	Read the User Manual before use.
	To protect the environment, dispose of the device at appropriate local collection points of electrical and electronic products. Do not dispose of it in household waste.
	Medical device needs to be protected from moisture.
	Temperature limits to which the medical device can be safely exposed (transport / storage conditions -10 to +45°C and operating conditions +16 to +40°C).
	The range of humidity to which the medical device can be safely exposed (transport / storage and operating conditions 10%-85%).
	The limits of atmospheric pressure to which the medical device can be safely exposed (transport / storage and operating conditions 860hPa – 1060 hPa).
	Single use only

6. Warranty and service

The Higo® device is delivered with the finest service and warranty policy of 24 months from the date of purchase. The warranty shall only be valid if accessories approved by HigoSense are used and the device is used as described in the User Manual and according to the intended use. Do not use the device if it has been damaged. Any actions connected with repairs or replacements may only be performed by personnel of an authorized service of HigoSense. The manufacturer is not liable for any damages caused by non-compliance with the user manual.

If you have any problem or issue with your device, firstly see Section 4, Troubleshooting, and if you do not find a solution, please contact your healthcare provider or the manufacturer support via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com. When sending the device for service, it must be delivered in the original manufacturer's packaging (including paper box).

7. Adverse events and medical incidents

Any adverse event or serious incident that occurs in relation to the Higo® device should be reported to the manufacturer via the contact form on the webpage www.higosense.com or by sending an email to support@higosense.com and to the competent authority of the country in which the user is established.

Erratum no. 2 to the HIGO® User Manual v.6.0 of June 1, 2021. and v.7.0 of February 24, 2022.

NOTES in section 3.1.1 and 3.2.1 including Wi-Fi specification:

NOTE: A Wi-Fi router connected to the Internet is required to set up the mobile application. **Your Wi-Fi router must be on a 2.4 GHz or 5 GHz frequencies.** For more details about Wi-Fi router compatibility, refer to Section 5.1.1.

[...]

3. How to use the Higo®

3.3 Medical exams

[...]

A. Temperature measurement

Keep the thermometer module clean and protect it from dust or other small contaminants for high accuracy and repeatability of measurements. In the event of contamination inside the aluminum cone or the temperature sensor ^(4.1), **follow the cleaning instructions in section 3.5.1c.** Do not use body or dish soap or detergents. [...]

C. Left and right ear exam

[...] Step 2: Apply the correct speculum size according to the patient's age, by inserting the speculum onto the otoscope tip and twisting it until it is locked into place.

Choose the large speculum marked with size "4.0 mm" ^(2.1) for patients above 3 years of age (>3).

Choose the small speculum marked with size "2.5 mm" ^(2.2) for patients **until the age of 3** (≤3).

[...]

D. Throat exam

[...] Step 2: Connect the correct tongue depressing tip according to the patient's age, by placing the tip into the throat module slot until it clicks and is locked into place automatically.

Choose the large tongue depressing tip marked with the letter "L" ^(5.1) for patients above 5 years old (>5).

Choose the small tongue depressing tip marked with the letter "S" ^(5.2) for patients **until the age of 5** (≤5). [...]

3.5 Cleaning and maintenance

3.5.1 Cleaning

c. The most delicate parts of the device to pay a special attention during cleaning are the following:

- the lenses located on the front of the otoscope module,
- the cone and sensor of the thermometer module

These parts must be clean and intact to ensure proper results can be obtained. **Wipe their surface very gently with a dry cloth. If contamination is still present wipe them with a cotton swab moistened with alcohol for at least 30 seconds.** Hold the modules upside down to prevent excess moisture from entering the lenses and sensor area. After the alcohol has dried out (minimum after 2 minutes), you can perform an exam.

4.Troubleshooting

Situation	Solution
I cannot add the Wi-Fi network	Make sure the Wi-Fi name and password you are entering match and are correct. If the name and password are correct, turn it off and on and retry the process of adding Wi-Fi. If the issue persists, please contact your healthcare provider or HigoSense support.
The measured temperature is inaccurate or permanently out of range due to small contaminants present on the temperature sensor ^(4.1) .	Follow the cleaning instructions in section 3.5.1c

5.1 Specification of the device

Wireless connection frequency band	2.4 GHz and 5 GHz (see Section 5.1.1)
------------------------------------	--

5.1.1 Wi-Fi network specification

The following Wi-Fi network specifications are supported by the Higo® device:

1. **IEEE 802.11n standard is recommended for optimal performance.**
IEEE 802.11g standard is also supported.
Please note, that both bands 2.4 GHz and 5GHz are supported.
This device within the 5150-5250 MHz band is restricted to indoor use only.
2. **Transmitter technical characteristics:**
Frequency band: 2400 MHz – 2483.5 MHz, 5150-5250 MHz, 5250-5350 MHz and 5470-5725 MHz
Maximum output power:
20 dBm (100 mW) **for 2400 MHz – 2483.5 MHz**
23 dBm for 5150-5250 MHz
23 dBm for 5250-5350 MHz
30 dBm for 5470-5725 MHz
3. **Antenna information:**
Type: internal, dual band 2.4 and 5 GHz
Gain: 0.712 dBi **for 2.4 GHz** and 1.25 dBi **for 5 GHz.**

5.2 Declaration of Conformity

The safety, reliability and performance of this device can only be assured under the following conditions:

- a) The device has been used for its indicated use and according to the operating instructions in this User Manual.
- b) All fittings, extensions, readjustments, changes, or repairs have been carried out by HigoSense.

HigoSense declares that the Higo® conforms with the essential requirements and other relevant provisions of the **Council Directive 93/42/EEC**, the **Directive 2011/65/EU of the European Parliament and of the Council** and the **Directive 2014/53/UE of the European Parliament and of the Council**.

5.3 Description of label and packaging marking

Symbol	Description
	Manufacturer
	CE marking that indicates that a product complies with the Directive 93/42/EEC . Certified product. Notified body number: 2274

6. Guarantee and service

Higo® device is covered with Manufacturer's Guarantee binding for the period of 12-months from the date of the sales of the device. Detailed guarantee conditions can be found in the Guarantee Card enclosed with the Higo® device or the agreement. The manufacturer's guarantee is only granted to entities that are entrepreneurs (B2B) or legal entities.

For maintenance and repair service contact your **HEALTHCARE PROVIDER** or the manufacturer support via contact form on higosense.com or mail support@higosense.com.

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